



# Mobi Mesh

## Portable Nebulizer with Mesh Technology

APEX Mobi Mesh is an innovative, lightweight nebulizer combines the best of everything you want—small size, virtually silent, easy to use, and fine particle size that runs and operates anywhere! It is an ideal for use in children and it may “Put Easy Respiratory in your Hands Anytime”.

# Mobi Mesh

## Portable Nebulizer with Mesh Technology

### Features ▼



#### Mesh Technology

Advanced technology offering fine mist at  $<5\mu\text{m}$  (MMAD) with evenly distributed particle sizes.



#### Turbo Mode

Enable vaporization at turbo speed for in-depth cleaning.



#### Ideal for Children

Virtually silent operation and continuous vaporization at various angles.



#### Hygienic & Thoroughly Clean

Medical cup can be easily removed by 1-click & thoroughly cleaned with hot water.



#### Lightweight & Handy

Tubing-Free handy design that weigh only 138g, it is ideal for use anywhere.



#### Friendly to Medicine

Vaporize under a heat-less & water-free interface with less influence to the medicine.



#### Low Medication Residue

Vaporize until the last drop to avoid any wasted medicine.



#### Durable Mesh Design

Equipped with stainless steel mesh for maximum durability.




Hygienic medical cup with easy removal and resistant to hot water cleaning.



Tubing free design with silent operation and effective drug delivery.



Include 1 mouthpiece, 1 adult mask and 1 pediatric mask.

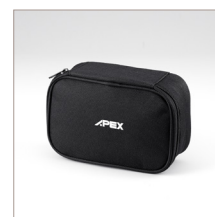

Include 1 carrying bag for easy storage and hand carry.

### Specification

Model Number	9R-001
Power Supply	1.5VX 2, AA size alkaline battery
Power Consumption	$< 2\text{W}$
Vibrating Frequency	103~123 kHz $\pm 1\text{K}$
Battery life	Up to 1.5 hrs. 2 weeks if daily use for 6 min treatment cycle
Flow Rate	0.35 $\pm 0.1$ mL/min
Particle Size	$< 5\mu\text{m}$
AUTO-OFF Time	20 min
Medical Cup Volume	8mL (max.)
Weight	138g (exclude battery)
Dimension	44.5mm (L) x 53.4mm (W) x 121.8mm (H)

#### Apex Medical Corp.

No. 9, Min Sheng St., Tu-Cheng,  
New Taipei City, 23679, Taiwan  
T +886 2 2268 3100  
F +886 2 2268 6525

marketing@apexmedicalcorp.com

www.apexmedicalcorp.com

Print-2014/05

All rights reserved Specifications may change without notice.